

X-ray photoelectron spectroscopy - An introduction


Spyros Diplas

spyros.diplas@sintef.no

spyros.diplas@smn.uio.no

*SINTEF Materials & Chemistry, Materials Physics -Oslo
&
Centre of Materials Science and Nanotechnology, Department of Chemistry, UiO*

Material Characterisation Methods

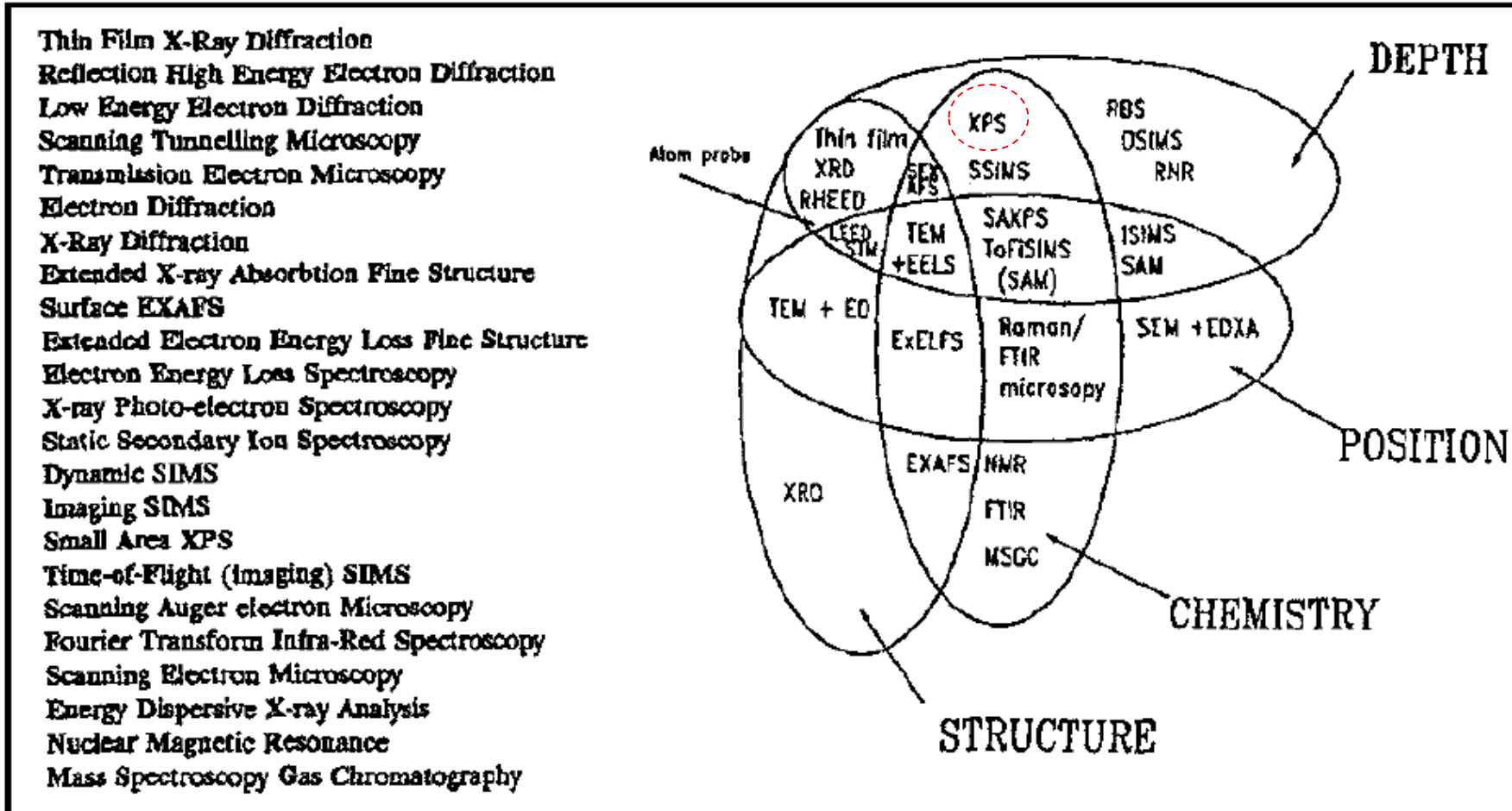

What is surface?


- What happens at surfaces is extremely important in a vast range of applications from environmental corrosion to medical implants.
- A surface is really the interface between different phases (solid, liquid or gas).
- We can think of the surface as the top layer of atoms but in reality the state of this layer is very much influenced by the 2 – 10 atomic layers below it (~0.5 – 3 nm).
- Surface modification treatments are often in the range of 10 – 100 nm thick. >100 nm can be thought of as the bulk.
- Surface analysis encompasses techniques which probe the properties in all these ranges.

*God made solids, but surfaces were the work of the devil
-----Wolfgang Pauli*


Surface Analysis - Techniques Available

- Properties and reactivity of the surface will depend on:
 - bonding geometry of molecules to the surface
 - physical topography
 - chemical composition
 - chemical structure
 - atomic structure
 - electronic state

No one technique can provide all these pieces of information. However, to solve a specific problem it is seldom necessary to use every technique available.

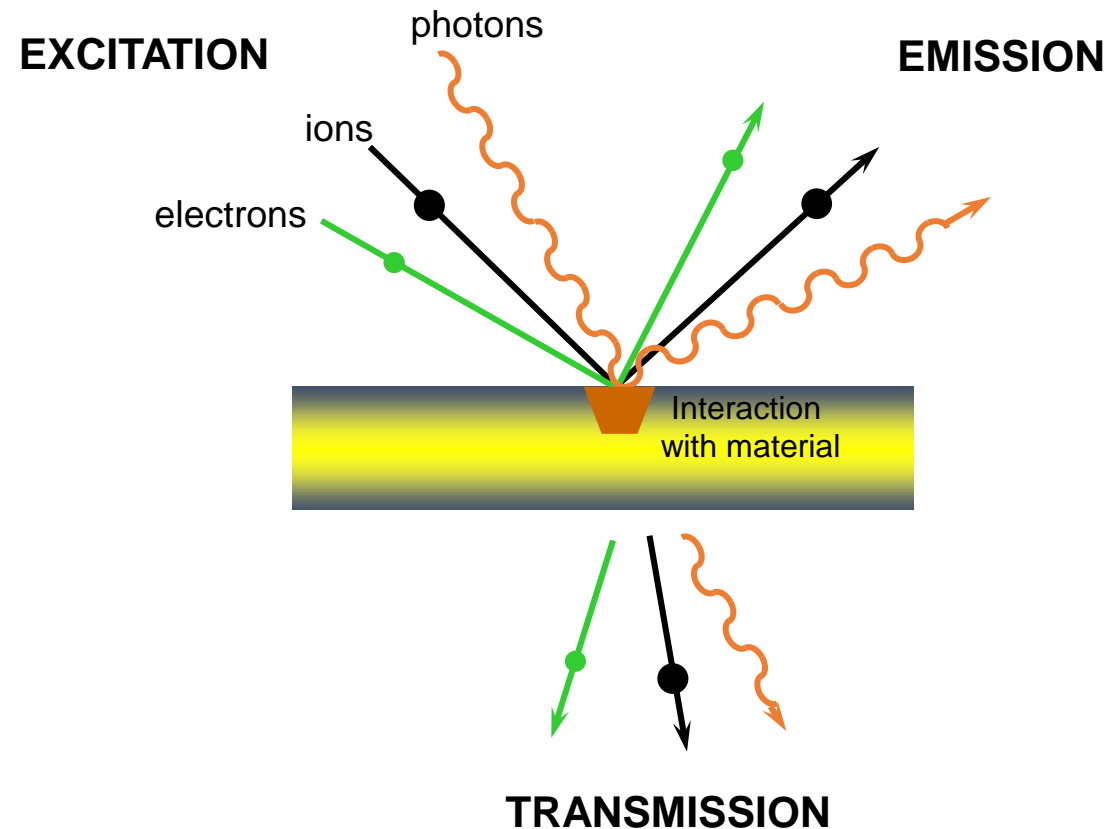


XPS-Basic Principle

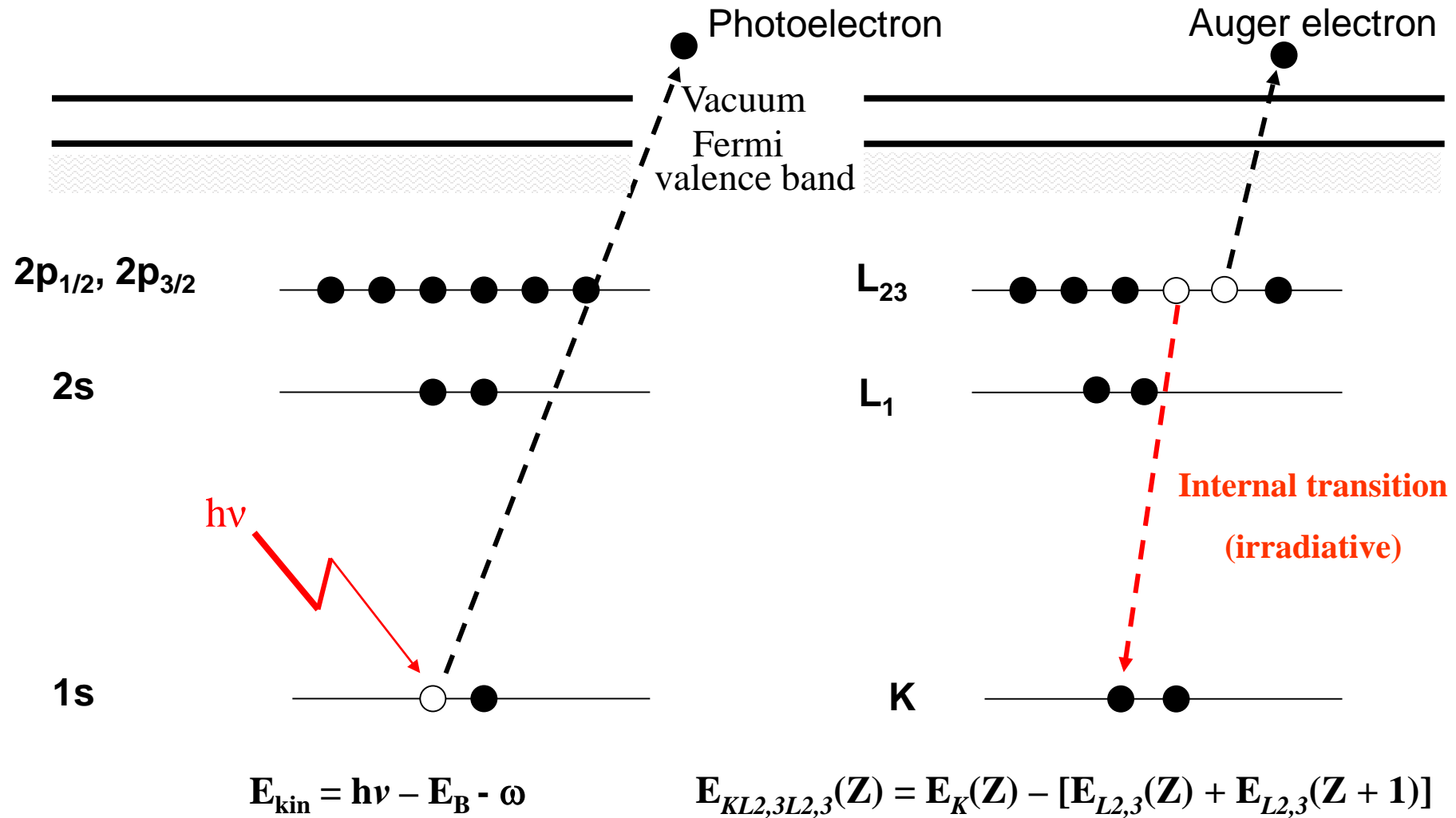


Excitation

De-excitation

Auger electron vs x-ray emission yield

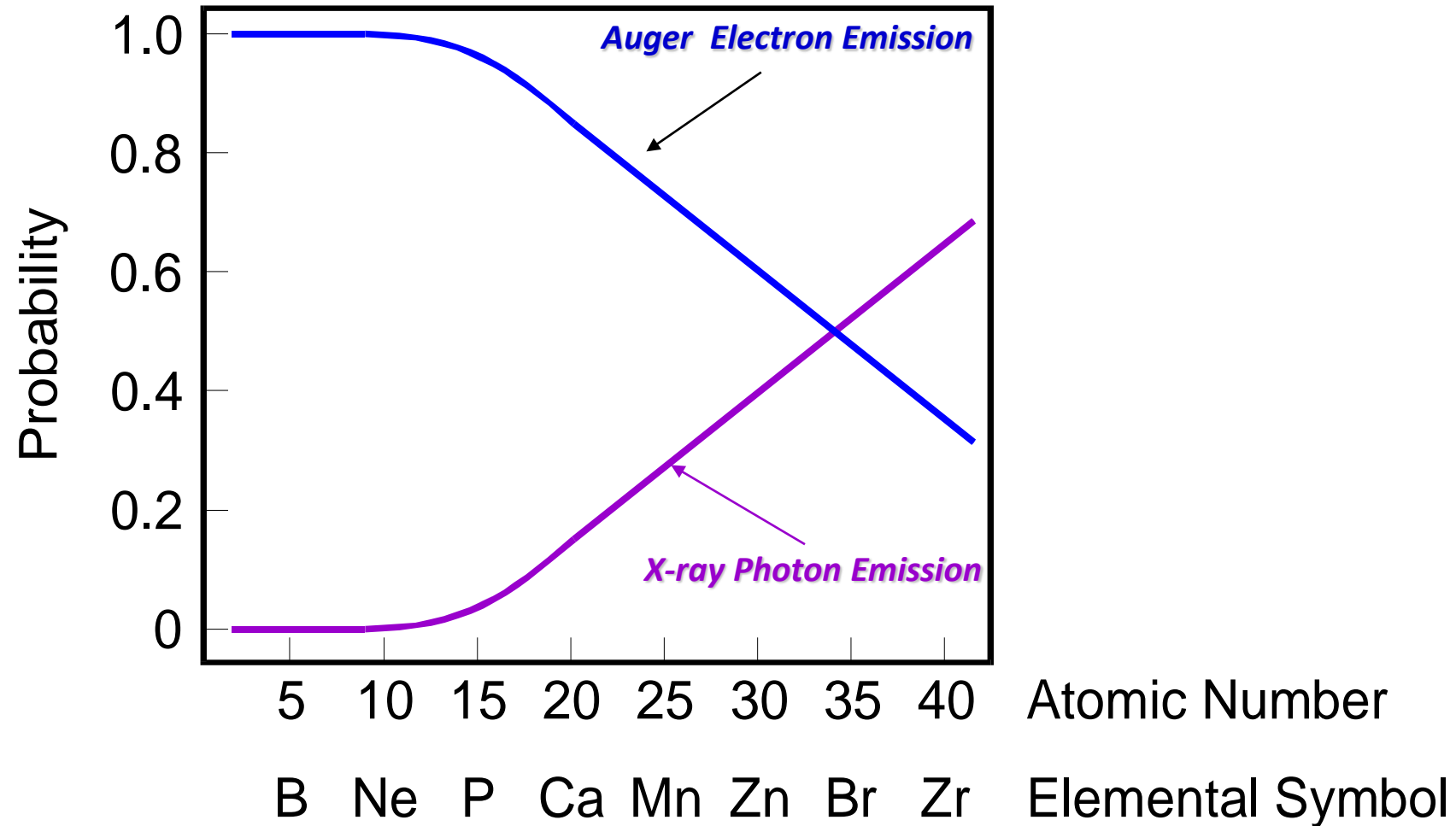

Schematic of an XPS spectrometer

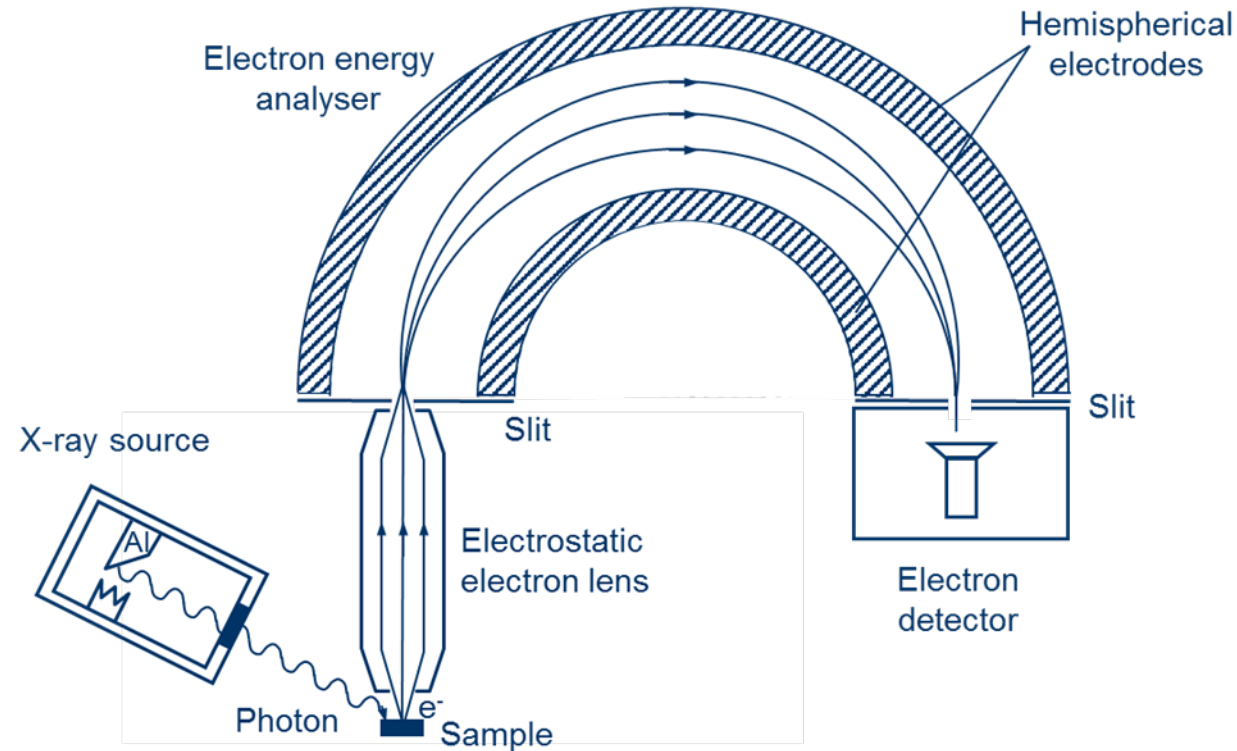


Number of emitted electrons measured as function of their kinetic energy

Examples of XPS spectrometers




Instrument: Kratos Axis Ultra^{DLD}

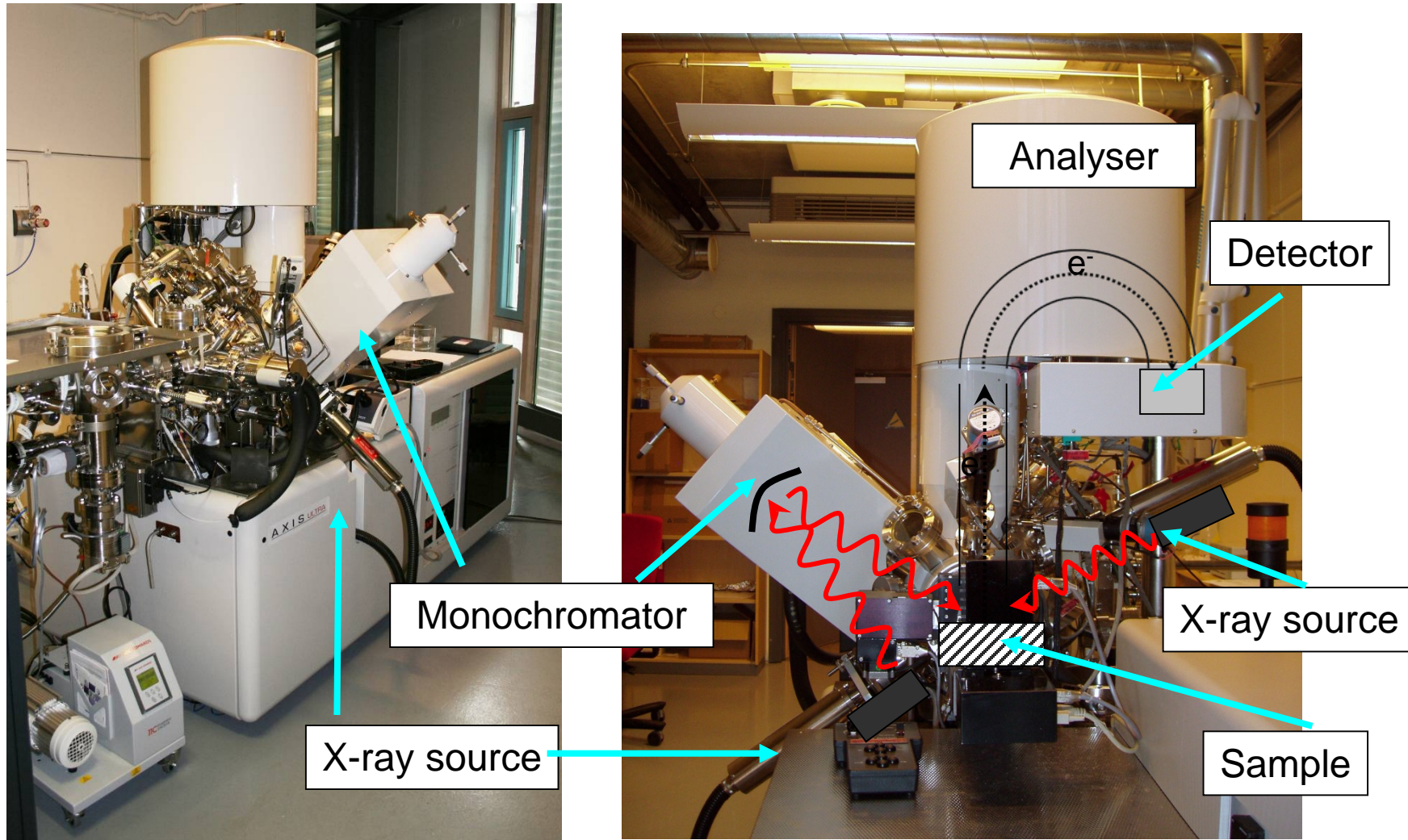

Instrument: Theta Probe

- Spectroscopy
 - **Source-defined** small area XPS
 - 15 μm to 400 μm
- Snapshot spectrum acquisition
 - Up to 112 channels
 - Faster serial mapping
 - Faster profiling
- **Unique parallel ARXPS** with up to 96 channels
- **Large samples** (70 mm x 70 mm x 25 mm)
- Sputter profiles
- Mapping possible up to full size of sample holder
- **ISS** included

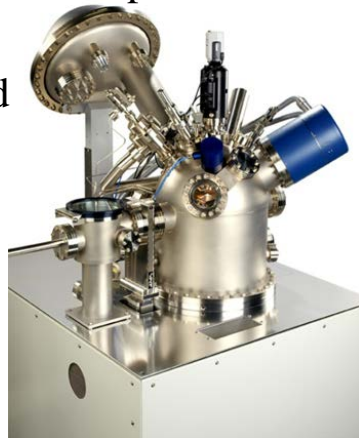

Target applications

- Thickness measurements
- Surface modification, plasma & chemical
- Self assembly
- Nanotechnology
- Ultra thin film technologies
- Shallow interfaces

Sample requirements


- Has to withstand high vacuum ($\leq 10^{-7}$ Torr).
- Has to withstand irradiation by X-rays
- Sample surface must be clean!
- Reasonably sized.

XPS Depth of Analysis

The probability that a photoelectron will escape from the sample without losing energy is regulated by the Beer-Lambert law:

$$e^{-\frac{z}{\lambda_e(E) \cos\theta}}$$

Where λ_e is the photoelectron inelastic mean free path

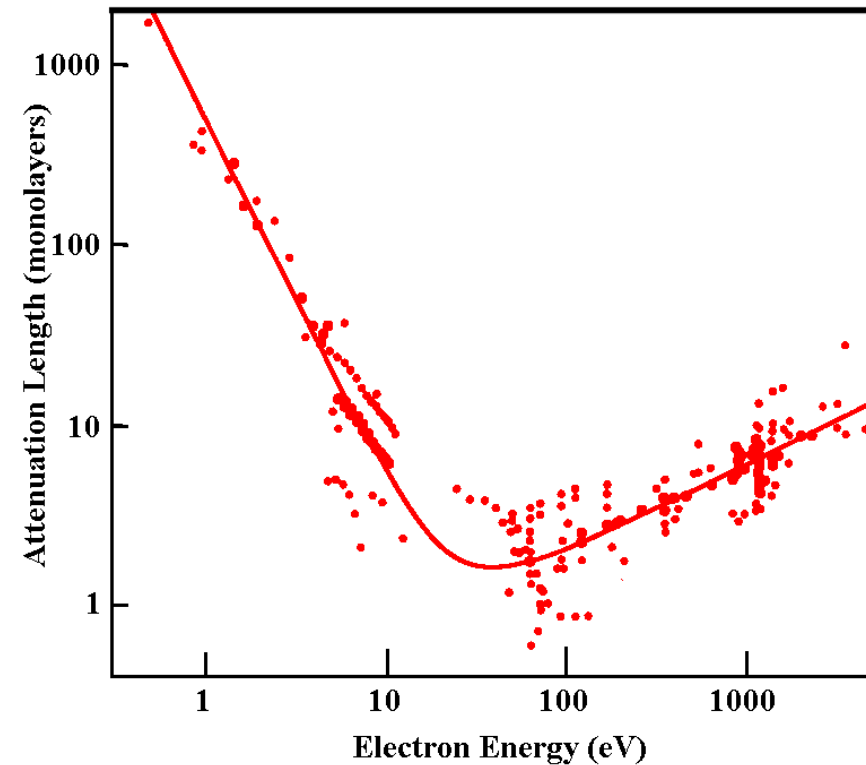

Attenuation length (λ) ≈ 0.9 IMFP

IMFP: The average distance an electron with a given energy travels between successive inelastic collisions

Features of the XPS spectrum


■ Primary structure

- Core level photoelectron peaks (atom excitation)
- Valence band spectra
- CCC, CCV, CVV Auger peaks (atom de-excitation)

■ Secondary structure

- X-ray satellites and ghosts
- Shake up and shake off satellites
- Plasmon loss features
- Background (slope)

XPS spectrum ITO

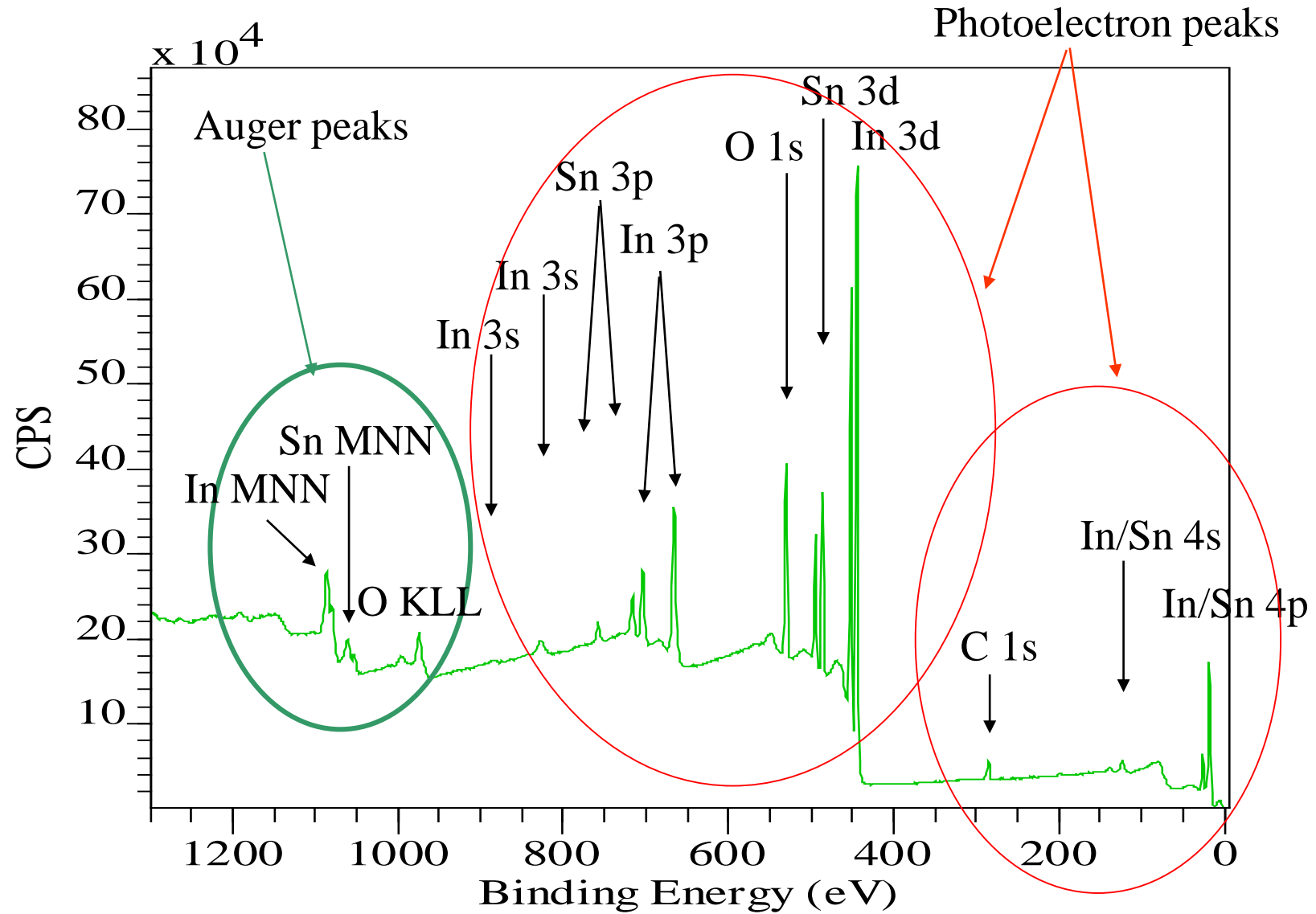

Peak width (ΔE)

$$\Delta E = (\Delta E_n^2 + \Delta E_p^2 + \Delta E_a^2)^{1/2}$$

Natural width X-ray source contribution Analyser contribution

- Gaussian broadening:

-Instrumental:

There is no perfectly resolving spectrometer nor a perfectly monochromatic X-ray source.

-Sample

For semiconductor surfaces in particular, variations in the defect density across the surface will lead to variations in the band bending and, thus, the work function will vary from point to point. This variation in surface potential produces a broadening of the XPS peaks.

-Excitation process such as the shake-up/shake-off processes or vibrational broadening.

- Lorentzian broadening.

The core-hole that the incident photon creates has a particular lifetime (τ) which is dependent on how quickly the hole is filled by an electron from another shell. From **Heisenberg's uncertainty principle**, the finite lifetime will produce a broadening of the peak.

$$\Gamma = \hbar/\tau$$


Intrinsic width of the same energy level should increase with increasing atomic number

Chemical shift

$$\Delta E_{(i)} = k\Delta q + \Delta V_M - \Delta R$$




Initial state contribution

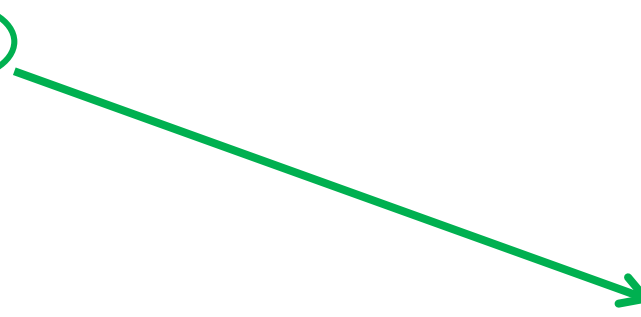

final state contribution

- Δq : changes in valence charge
- ΔV_M : Coulomb interaction between the photoelectron (i) and the surrounding charged atoms.
- ΔR : relaxation energy change arising from the response of the atomic environment (local electronic structure) to the screening of the core hole

Chemical shift - Growth of ITO on p c-Si

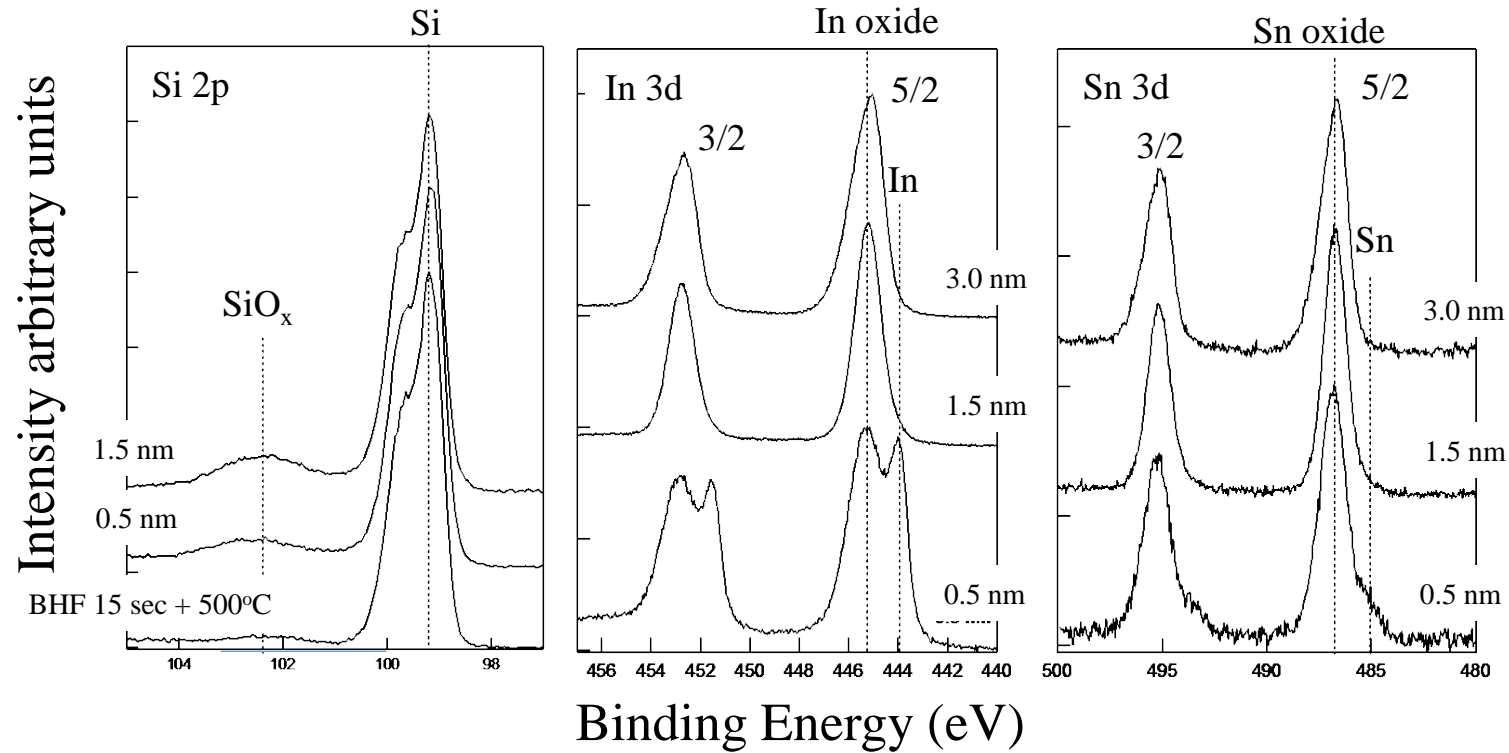



Quantification

- Unlike AES, SIMS, EDX, WDX there are little in the way of matrix effects to worry about in XPS. We can use either theoretical or empirical cross sections, corrected for transmission function of the analyser. In principle the following equation can be used:


$$I = J \rho \sigma K \lambda$$

- **I** is the electron intensity
 - **J** is the photon flux,
 - **ρ** is the concentration of the atom or ion in the solid,
 - **σ** is the cross-section for photoelectron production (which depends on the element and energy being considered),
 - **K** is a term which covers instrumental factors,
 - **λ** is the electron attenuation length.
-
- In practice atomic sensitivity factors (**F**) are often used:
 - $[A] \text{ atomic \%} = \{(IA/FA)/\Sigma(I/F)\}$
 - Various compilations are available.


Chemical shift

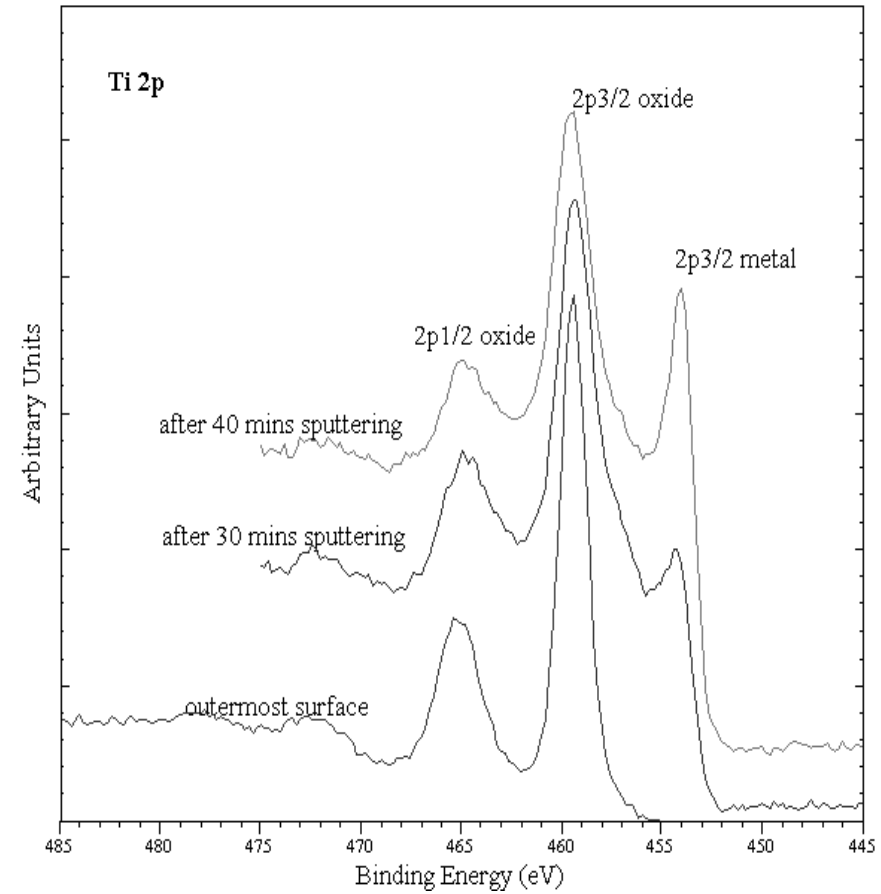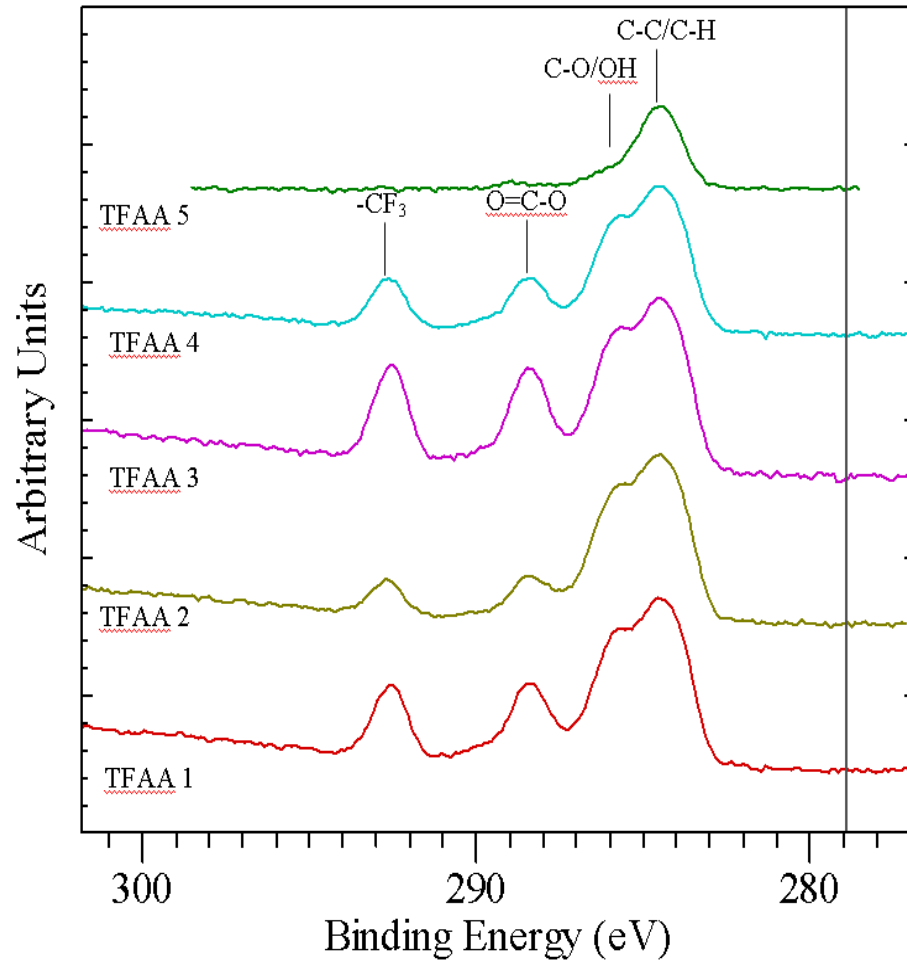

Spin-Orbit Coupling/Splitting

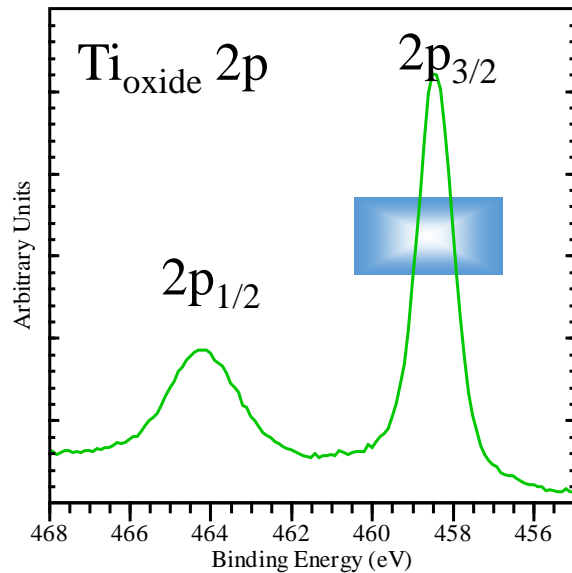

The intensity of the peaks is given by the degeneracy

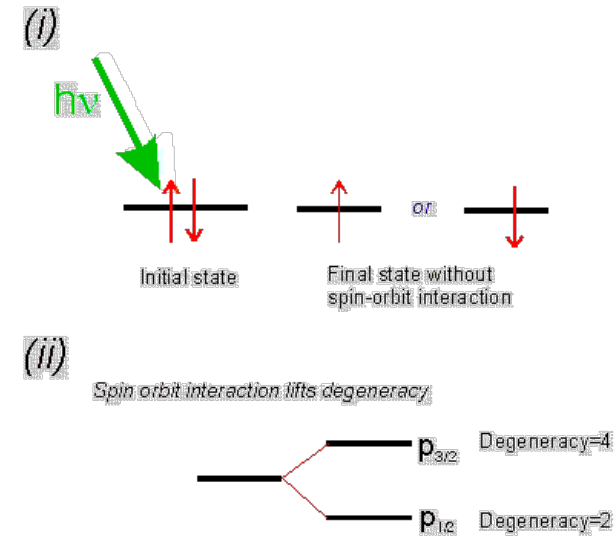
$$g_J = 2j+1$$


- **Spin-orbit coupling/ splitting:** final state effect for orbitals with orbital angular momentum $l > 0$. A **magnetic interaction** between an electron's spin and its orbital angular momentum.
- Example Ti. Upon photoemission an electron from the p orbital is removed - remaining electron can adopt one of two configurations: a *spin-up* ($s=+1/2$) or *spin-down* ($s=-1/2$) state. If no spin-orbit interaction these two states would have equal energy (degenerated states).
- spin-orbit coupling **lifts the degeneracy**
- To realise that we need to consider the quantum number, j , the *total angular momentum quantum number*.
- $j=l + s$ where s is the spin quantum number ($\pm 1/2$). For a p orbital $j=1/2$ or $3/2$. Thus the final state of the system may be either $p_{1/2}$ or $p_{3/2}$ and this gives rise to a splitting of the core-level into a doublet as shown in the figure above.
- **Spin-orbit coupling** is described for light elements by the Russell-Saunders (LS) coupling approximation and by the $j-j$ coupling approximation for heavier elements

Plasmons





- They describe the interaction (inelastic scattering) of the PE with the plasma oscillation of the outer shell (valence band) electrons
- Plasmons in their quantum mechanical description are pseudoparticles with energy $E_p = \hbar\omega$
- $\omega = (\mathbf{n}e^2/\epsilon_0 m)^{1/2}/2\pi$ \mathbf{n} = valence electron density,
 e , m electron charge and mass
 ϵ_0 = dielectric constant of vacuum

Pure elements

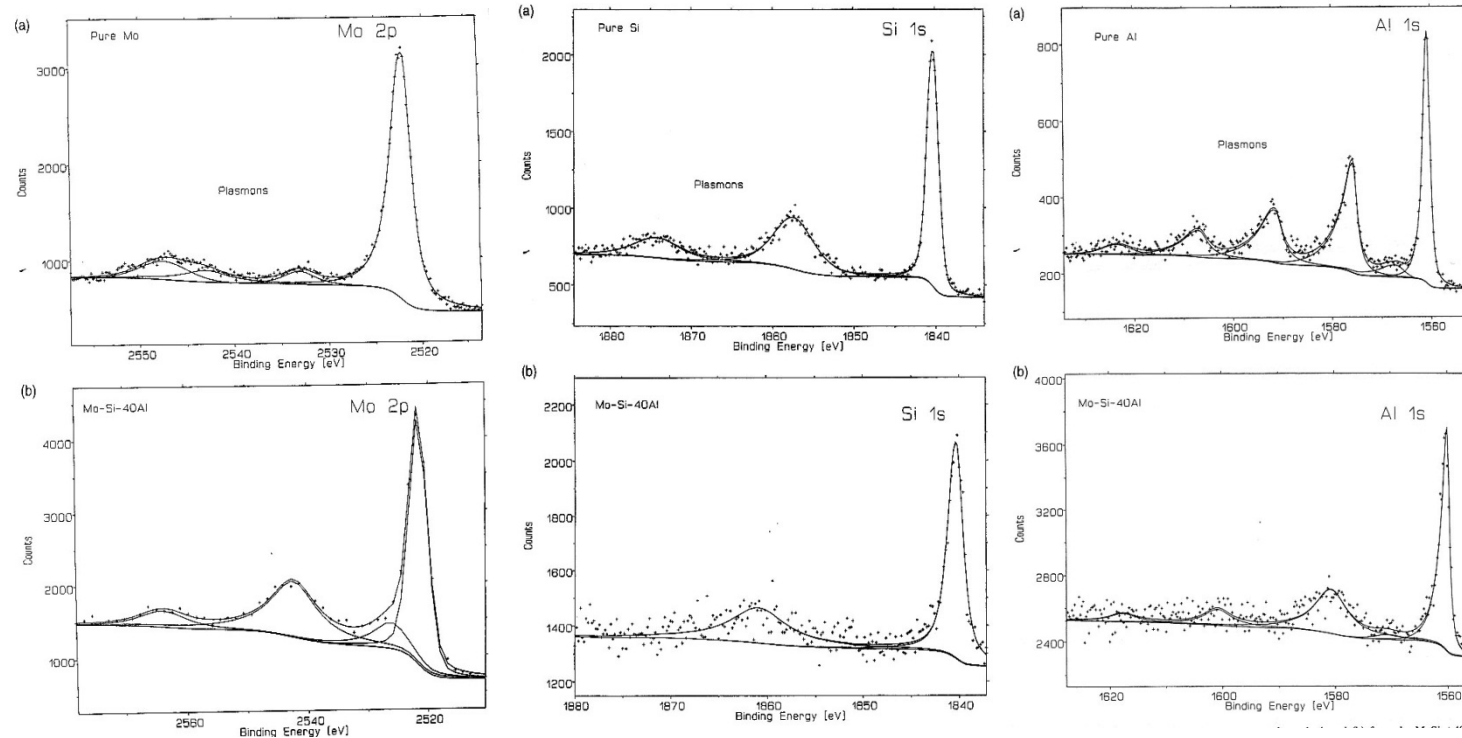


Mo-Si-Al
Compound

Peak asymmetry

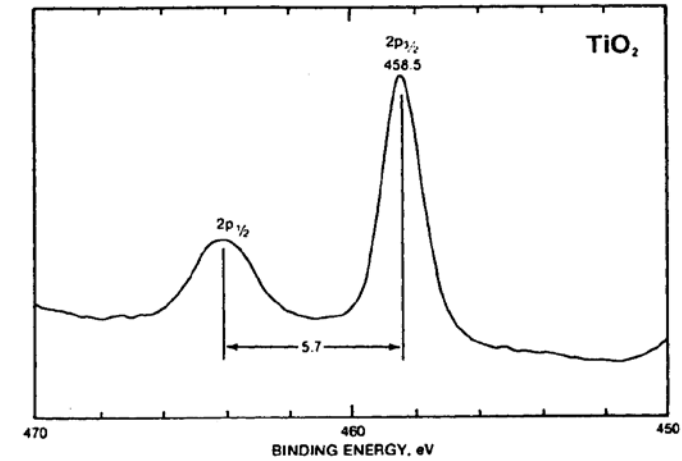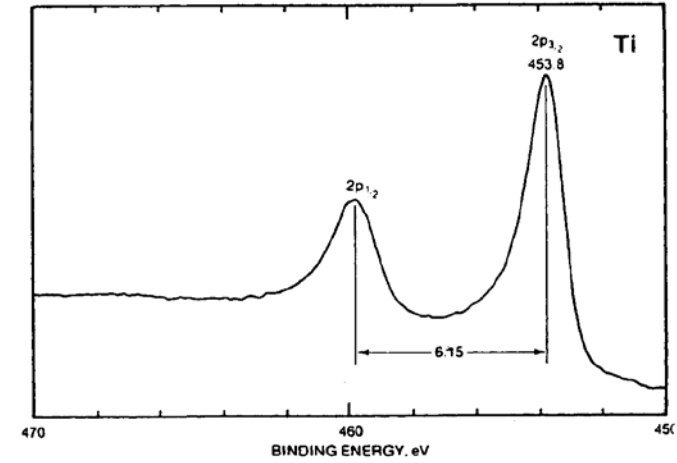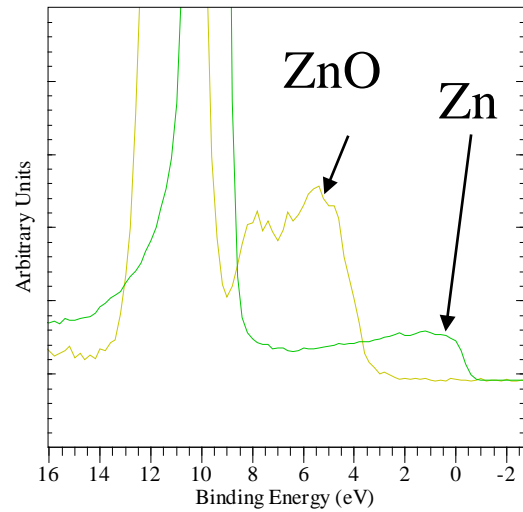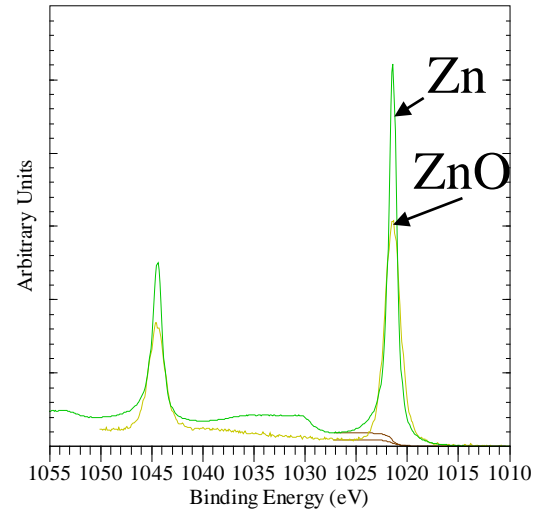


Peak asymmetry in metals caused by small energy electron-hole excitations near E_F of metal



Depth profile with ion sputtering

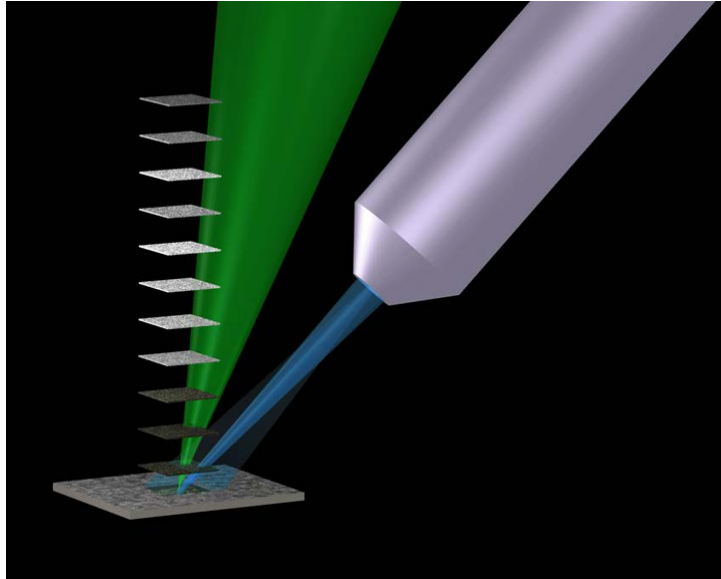

- Use of an ion gun to erode the sample surface and re-analyse
- Enables layered structures to be investigated
- Investigations of interfaces
- Depth resolution improved by:
 - Low beam energies
 - Small ion beam sizes
 - Sample rotation

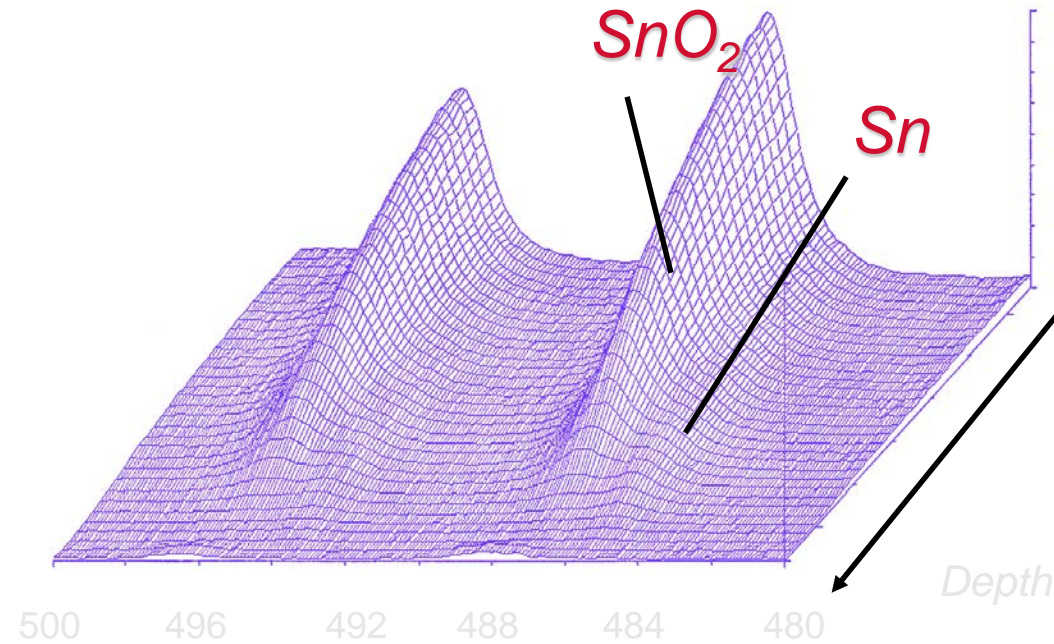


Angle Resolved XPS (ARXPS) for non-destructive depth profile

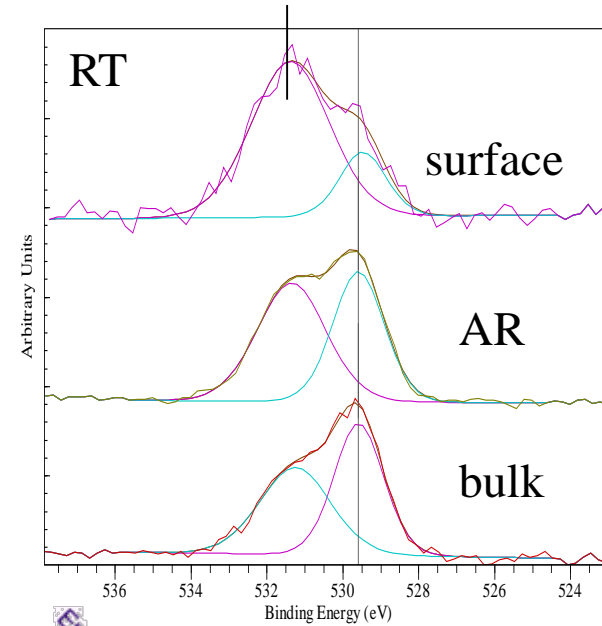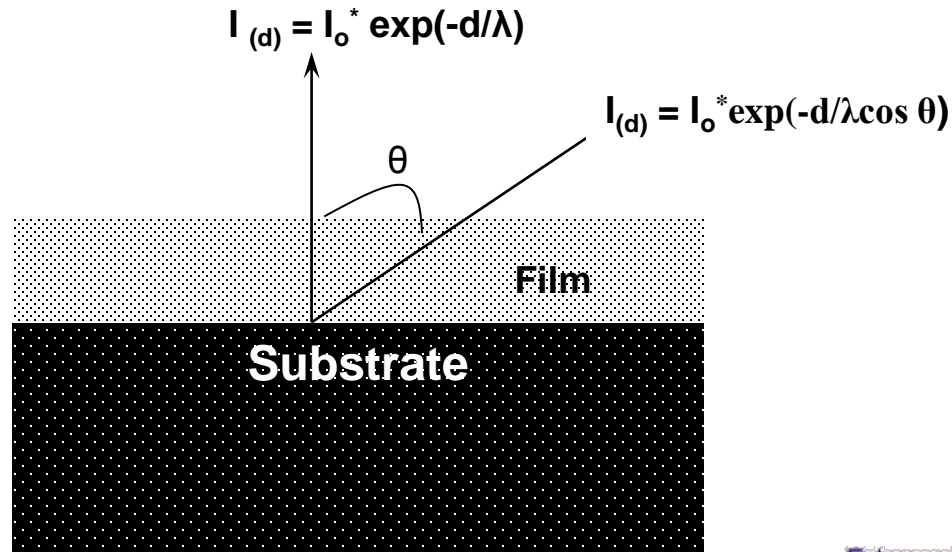

λ =attenuation length ($\lambda \approx 0.9$ IMFP)

$$\lambda = 538 \alpha_A / E_A^2 + 0.41 \alpha_A (\alpha_A E_A)^{0.5}$$

(α_A^3 volume of atom, E_A electron energy)

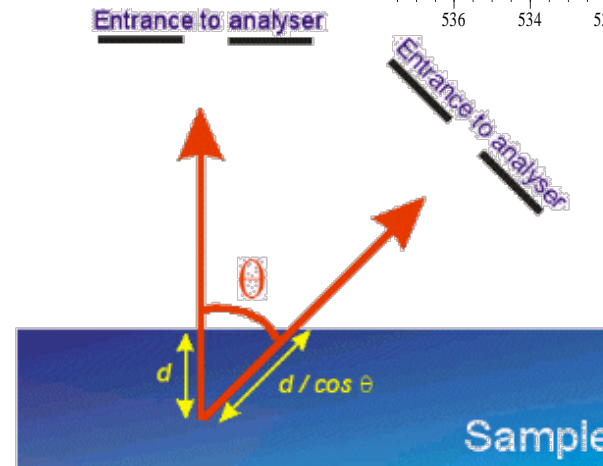

XPS-Check list

- Depth of analysis ~ 5nm
- All elements except H and He
- Readily quantified (limit ca. 0.1 at%)
- All materials (vacuum compatible)
- Chemical/electronic state information
 - Identification of chemical states*
 - Reflection of electronic changes to the atomic potential*
- Compositional depth profiling by
 - ARXPS (ultra thin film <10 nm),
 - change of the excitation energy
 - choose of different spectral areas
 - sputtering
- Ultra thin film thickness measurement
- Analysis area mm² to 10 micrometres


Interfacial studies of Al₂O₃ deposited on 4H-SiC(0001)

Avice, Diplas, Thøgersen, Christensen, Grossner, Svensson, Nilsen, Fjellvåg, Watts
Appl. Physics Letters, 2007;91, 52907, *Surface & Interface Analysis*, 2008;40,822


- $d = \lambda_{Si} \cos\theta \ln(1+R/R^\infty)$
- d: SiO_x film thickness
- λ_{Si} : inelastic mean free path for Si,
- θ : the angle of emission,
- R: the Si 2p intensity ratios I_{SiO_x}/I_{SiC},
- R[∞] the Si 2p intensity ratios I_{SiO_x}[∞]/I_{SiC}[∞] where I[∞] is the intensity from an infinitely thick substrate.
- $R^\infty = (\sigma_{Si, SiO_2} \cdot \lambda_{Si, SiO_2}) / (\sigma_{Si, Si} \cdot \lambda_{Si, Si})$
- where σ_{Si, SiO_2} and λ_{Si, SiO_2} are the number of Si atoms per SiO₂ unit volume and the inelastic mean free path respectively
- The $\sigma_{Si, SiO_2} / \sigma_{Si, Si}$ ratio is given by
- $\sigma_{Si, SiO_2} / \sigma_{Si, Si} = (D_{SiO_2} \cdot F_{Si}) / D_{Si} \cdot F_{SiO_2}$
- where D is the density of the material and F the formula weight.
- For the calculations we also assumed that the Si 2p photoelectrons from both SiC and Si oxide film will be attenuated by the same amount as they travel through the Al₂O₃ film therefore, their intensity ratio will reflect the attenuation of the Si 2p electrons coming from the SiC through the Si oxide film.

From XPS

d= 1nm at RT, d=3nm at 1273 K

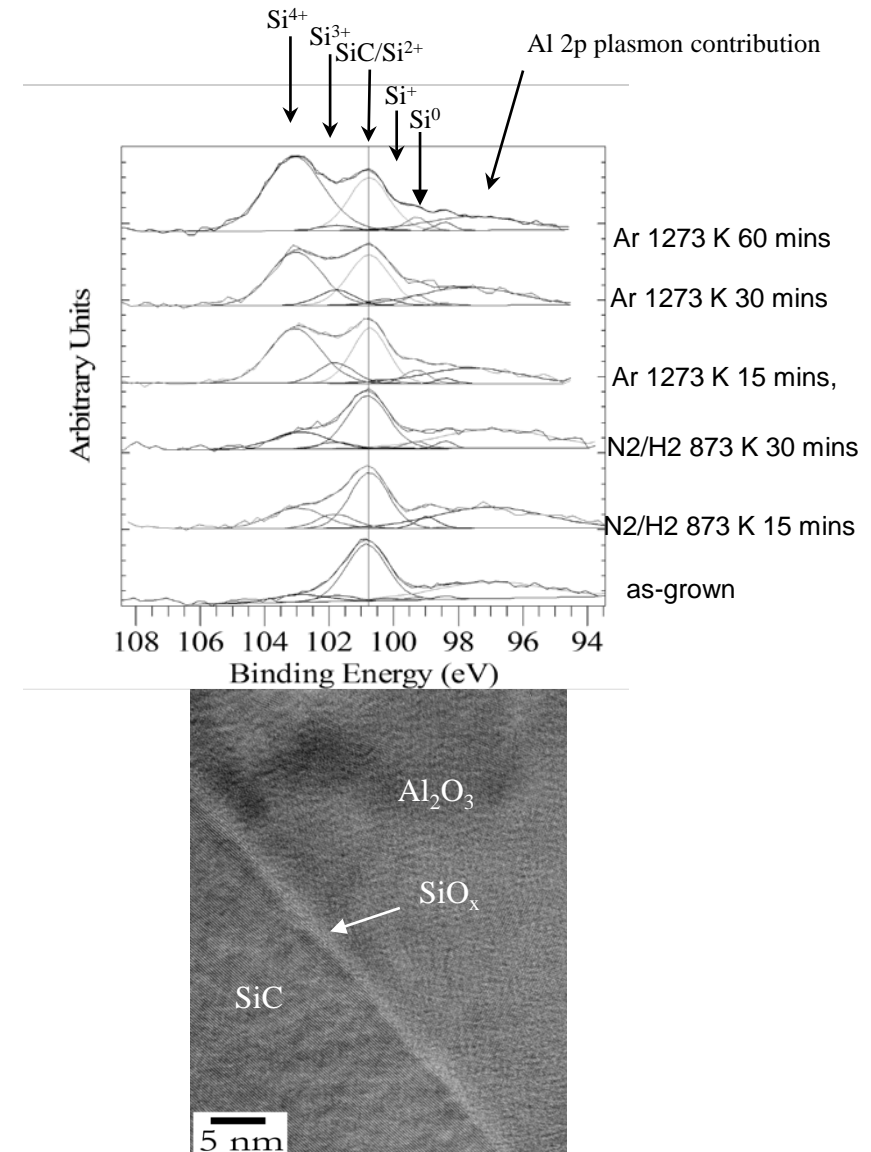

CIGS solar cell

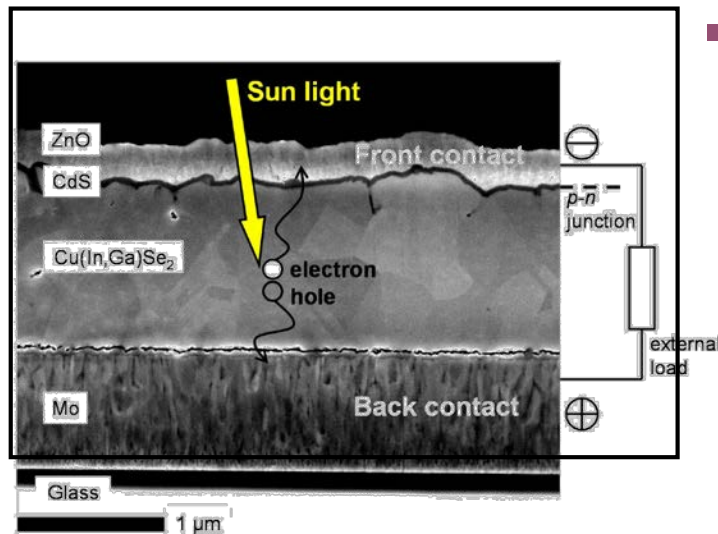

³SEM of a Cu(In,Ga)Se₂ solar cell (cross-section) and its mode of operation

■ CIGS solar cell

- **Energy/environmental application**
 - Solar cells based on Cu(In, Ga)Se₂ (CIGS)
 - Thin-film stack on glass
 - Mo and Zn oxide layer form electrical contacts
 - *p*-type CIGS film (sunlight absorber) and *n*-type CdS film form *p-n* junction
 - Excellent efficiency
 - Low cost compared to thicker silicon-based solar cells
- **Practical problem**
 - Controlling film composition and interfacial chemistry between layers (affects electrical properties)
- **XPS solution**
 - XPS sputter depth profiling
 - *Elemental and composition information as a function of depth*
 - *Identify chemical gradients within layers*
 - *Investigate chemistry at layer interfaces*


Acknowledgement: Thermo Electron Corporation

CIGS solar cell


■ CIGS solar cell

- Depth profile of CIGS film stack

- Demonstrates standardless quantification of XPS
 - Excellent quantification agreement between XPS and Rutherford BackScattering (RBS)
 - Both techniques show cross-over of In and Ga close to 1.6 μm depth
 - XPS tool is able to analyze product solar cell device

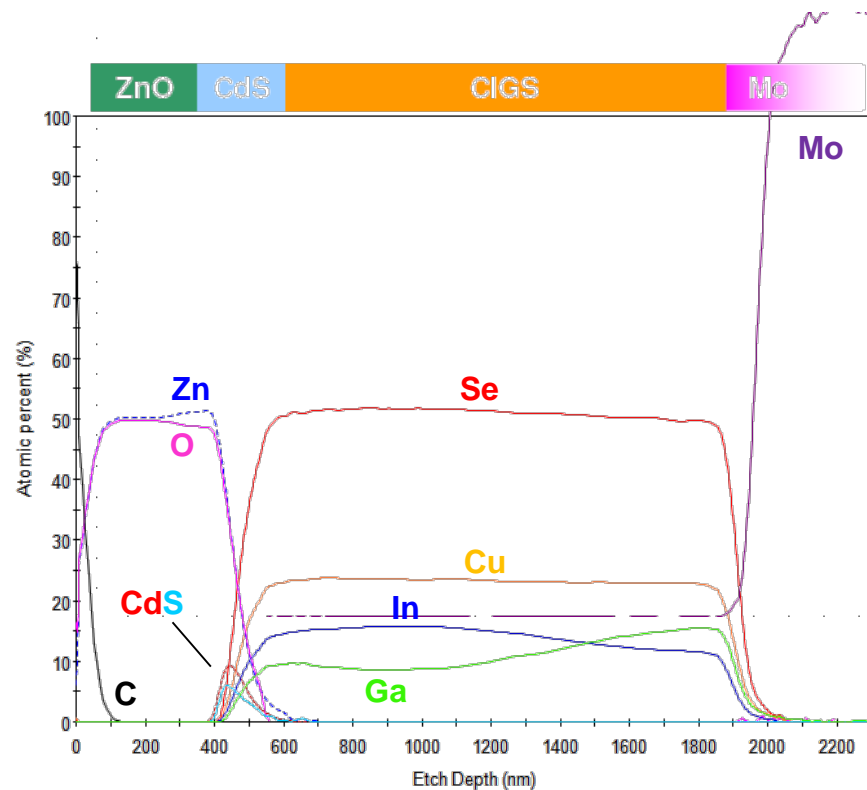

Sputter depth profile of CIGS film stack

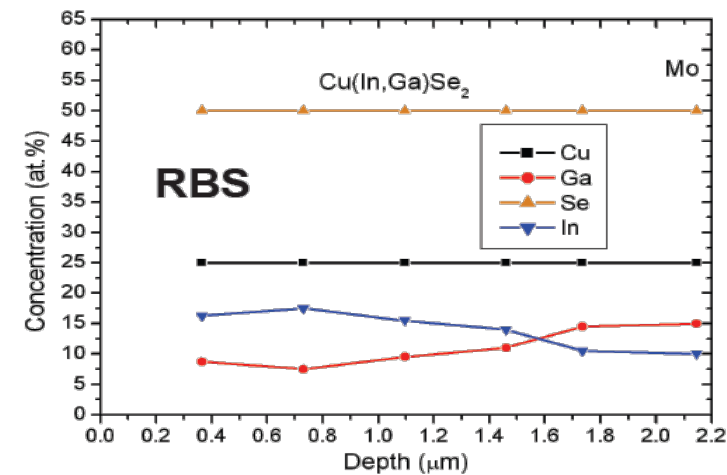

Rutherford BackScatter profile of CIGS film stack

Acknowledgement: Thermo Electron Corporation

Elemental distribution and oxygen deficiency of magnetron sputtered ITO films

A. Thøgersen, M.Rein, E. Monakhov, J. Mayandi, S. Diplac

JOURNAL OF APPLIED PHYSICS 109, 113532 (2011)

